

AESKULISA dsDNA-G

Ref 3142





Product Ref.	3142
Product Desc.	dsDNA-G
Manual Rev. No.	003 : 2013-10-10

Instruction Manual

Table of Contents

1	Intended Use	1
2	Clinical Application and Principle of the Assay	1
3	Kit Contents	2
4	Storage and Shelf Life	2
5	Precautions of Use.....	3
6	Sample Collection, Handling and Storage.....	4
7	Assay Procedure	4
8	Quantitative and Qualitative Interpretation	7
9	Technical Data.....	8
10	Performance Data.....	8
11	Literature	9





Product Ref.	3142
Product Desc.	dsDNA-G
Manual Rev. No.	003 : 2013-10-10

1 Intended Use

AESKULISA dsDNA-G is a solid phase enzyme immunoassay with human recombinant double-stranded DNA (dsDNA) for the quantitative and qualitative detection of IgG antibodies against dsDNA in human serum.

Anti-dsDNA antibodies mainly recognize the phosphat units of the DNA, thus these autoantibodies also bind single stranded DNA (ssDNA). To ensure correct quantitation of anti-dsDNA antibodies the used antigen has been proven to be free of contamination with ssDNA.

The assay is a tool in the differential diagnosis of systemic lupus erythematosus (SLE).

2 Clinical Application and Principle of the Assay

Antibodies binding to DNA belong to the group of anti-nuclear antibodies (ANA) that have been observed in several autoimmune diseases. Antibodies reacting with native double-stranded (ds) DNA are regarded as being specific for systemic lupus erythematosus (SLE) and have been observed in approximately 50-80% of the patients.

Antibodies against dsDNA are found during active phases of SLE. The amount of the serum concentration is positively correlated with the severity of the disease. Thus, detection of these autoantibodies is important for the diagnosis and the clinical monitoring of SLE. Consequently it has been established as 1 of the 11 ACR-criteria for the diagnosis of SLE.

Most patients with SLE display IgG class antibodies against dsDNA. These autoantibodies are associated with lupus nephritis. Approximately 30% of the SLE patients develop IgA class anti-dsDNA antibodies, additionally. There have been suggestions that the presence of these IgA class anti-dsDNA antibodies may define a certain subset of SLE patients. Indeed studies demonstrated the association of this subclass with certain parameters of the disease activity, such as elevated erythrocyte sedimentation rate, or the consumption of complement component C3, as well as the clinical parameters of cutaneous vasculitis, acral necrosis and erythema. While no association was found for nephritis and arthritis.

IgM class anti-dsDNA antibodies were found in 52 % of the sera from patients with SLE. In contrast to IgG and IgA class autoantibodies, the subclass IgM antibodies do not correlate with disease activity. However, a highly significant negative correlation between IgM anti-dsDNA antibodies and lupus nephritis, including its laboratory parameters was demonstrated. Therefore IgM class anti-dsDNA antibodies may indicate a subset of lupus patients being protected against the risk of developing nephritis.

Principle of the test

Serum samples diluted 1:101 are incubated in the microplates coated with the specific antigen. Patient's antibodies, if present in the specimen, bind to the antigen. The unbound fraction is washed off in the following step. Afterwards anti-human immunoglobulins conjugated to horseradish peroxidase (conjugate) are incubated and react with the antigen-antibody complex of the samples in the microplates. Unbound conjugate is washed off in the following step. Addition of TMB-substrate generates an enzymatic colorimetric (blue) reaction, which is stopped by diluted acid (color changes to yellow). The intensity of color formation from the chromogen is a function of the amount of conjugate bound to the antigen-antibody complex and this is proportional to the initial concentration of the respective antibodies in the patient sample.



3 Kit Contents

TO BE RECONSTITUTED				
Item	Quantity	Cap color	Solution color	Description / Contents
Sample Buffer (5x)	1 x 20ml	White	Yellow	5 x concentrated Tris, sodium chloride (NaCl), bovine serum albumin (BSA), sodium azide < 0.1% (preservative)
Wash Buffer (50x)	1 X 20ml	White	Green	50 x concentrated Tris, NaCl, Tween 20, sodium azide < 0.1% (preservative)
READY TO USE				
Item	Quantity	Cap color	Solution color	Description / Contents
Negative Control	1 x 1.5ml	Green	Colorless	Human serum (diluted), bovine serum albumin (BSA), sodium azide < 0.1% (preservative)
Positive Control	1 x 1.5ml	Red	Yellow	Human serum (diluted), bovine serum albumin (BSA), sodium azide < 0.1% (preservative)
Cut-off Calibrator	1 x 1.5ml	Blue	Yellow	Human serum (diluted), bovine serum albumin (BSA), sodium azide < 0.1% (preservative)
Calibrators	6 x 1.5ml	White	Yellow *	Concentration of each calibrator: 0, 3, 10, 30, 100, 300 IU/ml. Human serum (diluted), bovine serum albumin (BSA), sodium azide < 0.1% (preservative)
Conjugate, IgG	1 x 15ml	Blue	Blue	Containing: Anti-human immunoglobulins conjugated to horseradish peroxidase, bovine serum albumin (BSA)
TMB Substrate	1 x 15ml	Black	Colorless	Stabilized tetramethylbenzidine and hydrogen peroxide (TMB/H ₂ O ₂)
Stop Solution	1 x 15ml	White	Colorless	1M Hydrochloric Acid
Microtiter plate	12 x 8 well strips	N/A	N/A	With breakaway microwells. Refer to paragraph 1 for coating.
* Color increasing with concentration				
MATERIALS REQUIRED, BUT NOT PROVIDED				
Microtiter plate reader 450 nm reading filter and recommended 620 nm reference filter (600-690 nm). Glass ware (cylinder 100-1000ml), test tubes for dilutions. Vortex mixer, precision pipettes (10, 100, 200, 500, 1000 µl) or adjustable multipipette (100-1000µl). Microplate washing device (300 µl repeating or multichannel pipette or automated system), adsorbent paper. Our tests are designed to be used with purified water according to the definition of the United States Pharmacopeia (USP 26 - NF 21) and the European Pharmacopeia (Eur.Ph. 4th ed.).				

4 Storage and Shelf Life

Store all reagents and the microplate at 2-8°C/35-46°F, in their original containers. Once prepared, reconstituted solutions are stable at 2-8°C/35-46°F for at least 1 month. Reagents and the microplate shall be used within the expiry date indicated on each component, only. Avoid intense exposure of TMB solution to light. Store microplates in designated foil, including the desiccant, and seal tightly.



Product Ref.	3142
Product Desc.	dsDNA-G
Manual Rev. No.	003 : 2013-10-10

5 Precautions of Use

5.1 Health hazard data

THIS PRODUCT IS FOR IN VITRO DIAGNOSTIC USE ONLY. Thus, only staff trained and specially advised in methods of in vitro diagnostics may perform the kit. Although this product is not considered particularly toxic or dangerous in conditions of the intended use, refer to the following for maximum safety:

Recommendations and precautions

This kit contains potentially hazardous components. Though kit reagents are not classified being irritant to eyes and skin we recommend to avoid contact with eyes and skin and wear disposable gloves.

WARNING ! Calibrators, Controls and Buffers contain sodium azide (NaN_3) as a preservative. NaN_3 may be toxic if ingested or adsorbed by skin or eyes. NaN_3 may react with lead and copper plumbing to form highly explosive metal azides. On disposal, flush with a large volume of water to prevent azide build-up. Please refer to decontamination procedures as outlined by CDC or other local/national guidelines.

Do not smoke, eat or drink when manipulating the kit. Do not pipette by mouth.

All human source material used for some reagents of this kit (controls, standards e.g.) has been tested by approved methods and found negative for HbsAg, Hepatitis C and HIV 1. However, no test can guarantee the absence of viral agents in such material completely. Thus handle kit controls, standards and patient samples as if capable of transmitting infectious diseases and according to national requirements.

The kit contains material of animal origin as stated in the table of contents, handle according to national requirements.

5.2 General directions for use

In case that the product information, including the labeling, is defective or incorrect please contact the manufacturer or the supplier of the test kit.

Do not mix or substitute Controls, Calibrators, Conjugates or microplates from different lot numbers. This may lead to variations in the results.

Allow all components to reach room temperature (20-32°C/68-89.6°F) before use, mix well and follow the recommended incubation scheme for an optimum performance of the test.

Incubation: We recommend test performance at 30°C/86°F for automated systems.

Never expose components to higher temperature than 37°C/ 98.6°F.

Always pipette substrate solution with brand new tips only. Protect this reagent from light. Never pipette conjugate with tips used with other reagents prior.

A definite clinical diagnosis should not be based on the results of the performed test only, but should be made by the physician after all clinical and laboratory findings have been evaluated. The diagnosis is to be verified using different diagnostic methods.



6 Sample Collection, Handling and Storage

Use preferentially freshly collected serum samples. Blood withdrawal must follow national requirements. Do not use icteric, lipemic, hemolysed or bacterially contaminated samples. Sera with particles should be cleared by low speed centrifugation (<1000 x g). Blood samples should be collected in clean, dry and empty tubes.

After separation, the serum samples should be used during the first 8h, respectively stored tightly closed at 2-8°C/35-46°F up to 48h, or frozen at -20°C/-4°F for longer periods

7 Assay Procedure

7.1 Preparations prior to starting

Dilute concentrated reagents:

Dilute the concentrated sample buffer 1:5 with distilled water (e.g. 20 ml plus 80 ml).

Dilute the concentrated wash buffer 1:50 with distilled water (e.g. 20 ml plus 980 ml).

To avoid mistakes we suggest to mark the cap of the different calibrators.

Samples:

Dilute serum samples 1:101 with sample buffer (1x)

e.g. 1000 µl sample buffer (1x) + 10 µl serum. Mix well !

Washing:

Prepare 20 ml of diluted wash buffer (1x) per 8 wells or 200 ml for 96 wells

e.g. 4 ml concentrate plus 196 ml distilled water.

Automated washing:

Consider excess volumes required for setting up the instrument and dead volume of robot pipette.

Manual washing:

Discard liquid from wells by inverting the plate. Knock the microwell frame with wells downside vigorously on clean adsorbent paper. Pipette 300 µl of diluted wash buffer into each well, wait for 20 seconds. Repeat the whole procedure twice again.

Microplates:

Calculate the number of wells required for the test. Remove unused wells from the frame, replace and store in the provided plastic bag, together with desiccant, seal tightly (2-8°C/35-46°F).

7.2 Pipetting Scheme

We suggest pipetting calibrators, controls and samples as follows:

For <i>QUANTITATIVE</i> interpretation					For <i>QUALITATIVE</i> interpretation				
	1	2	3	4...		1	2	3	4...
A	Cal A	Cal E	P1		A	NC	P2		
B	Cal A	Cal E	P1		B	NC	P2		
C	Cal B	Cal F	P2		C	CC	P3		
D	Cal B	Cal F	P2		D	CC	P3		
E	Cal C	PC	P3		E	PC	...		
F	Cal C	PC	P3		F	PC	...		
G	Cal D	NC	...		G	P1	...		
H	Cal D	NC	...		H	P1	...		

CalA: calibrator A

CalB: calibrator B

CalC: calibrator C

CalD: calibrator D

CalE: calibrator E

CalF: calibrator F

PC: positive control

NC: negative control




CC: cut-off calibrator

P1: patient 1

P2: patient 2

P3: patient 3

7.3 Test Steps

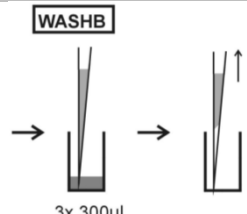Step	Description
1.	Ensure preparations from step 7.1 above have been carried out prior to pipetting.
2.	Use the following steps in accordance with quantitative/ qualitative interpretation results desired:
CONTROLS & SAMPLES	
3.	 <p>Pipette into the designated wells as described in chapter 7.2 above, 100 µl of either:</p> <ol style="list-style-type: none"> Calibrators (CAL.A to CAL.F) for <i>QUANTITATIVE</i> or Cut-off Calibrator (CC) for <i>QUALITATIVE</i> interp. <p>and 100 µl of each of the following:</p> <ul style="list-style-type: none"> Negative control (NC) and Positive control (PC), and Patients diluted serum (P1, P2...)
4.	 <p>Incubate for 30 minutes at 20-32°C/68-89.6°F.</p>
5.	 <p>Wash 3x with 300 µl washing buffer (diluted 1:50).</p>



CONJUGATE


6.

CONJ



Pipette 100 µl conjugate into each well.

7.

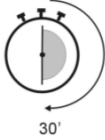

Incubate for 30 minutes at 20-32°C/68-89.6°F.

8.

WASHB



Wash 3x with 300 µl washing buffer (diluted 1:50).

SUBSTRATE

9.

SUB



Pipette 100 µl TMB substrate into each well.

10.




Incubate for 30 minutes at 20-32°C/68-89.6°F, protected from intense light.

STOP

11.

STOP

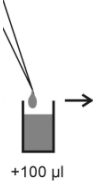

Pipette 100 µl stop solution into each well, using the same order as pipetting the substrate.

12.




Incubate 5 minutes minimum.

13.

Agitate plate carefully for 5 sec.

14.

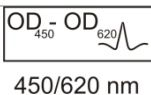

Read absorbance at 450 nm (recommended 450/620 nm) within 30 minutes.



8 Quantitative and Qualitative Interpretation

For **quantitative interpretation** establish the standard curve by plotting the **optical density (OD) of each calibrator (y-axis)** with respect to the corresponding concentration values in IU/ml (x-axis). For best results we recommend log/lin coordinates and 4-Parameter Fit. From the OD of each sample, read the corresponding antibody concentrations expressed in IU/ml.

Normal Range	Equivocal Range	Positive Results
< 12 IU/ml	12 - 18 IU/ml	>18 IU/ml

Example of a standard curve

Do NOT use this example for interpreting patient's result

Calibrators IgG	OD 450/620 nm	CV % (Variation)
0 IU/ml	0.036	2.9
3 IU/ml	0.176	2.3
10 IU/ml	0.314	2.9
30 IU/ml	0.618	2.9
100 IU/ml	1.312	0.1
300 IU/ml	2.076	0.7

Example of calculation

Patient	Replicate (OD)	Mean (OD)	Result (IU/ml)
P 01	0.799/0.744	0.772	40.3
P 02	1.404/1.393	1.399	119.5

Samples above the highest calibrator range should be reported as >Max. They should be diluted as appropriate and re-assayed. Samples below calibrator range should be reported as < Min.

For lot specific data, see enclosed quality control leaflet. Medical laboratories might perform an in-house quality control by using own controls and/or internal pooled sera, as foreseen by national regulations.

Each laboratory should establish its own normal range based upon its own techniques, controls, equipment and patient population according to their own established procedures.

In case that the values of the controls do not meet the criteria the test is invalid and has to be repeated.

The following technical issues should be verified: Expiration dates of (prepared) reagents, storage conditions, pipettes, devices, photometer, incubation conditions and washing methods.

If the items tested show aberrant values or any kind of deviation or that the validation criteria are not met without explicable cause please contact the manufacturer or the supplier of the test kit.

For **qualitative interpretation** read the optical density of the cut-off calibrator and the patient samples. Compare patient's OD with the OD of the cut-off calibrator. For qualitative interpretation we recommend to consider sera within a range of 20% around the cut-off value as equivocal. All samples with higher ODs are considered positive, samples with lower ODs are considered negative.

Negative:		OD patient	<	0.8 x OD cut-off	
Equivocal:	0.8 x	OD cut-off	≤	OD patient	≤ 1.2 x OD cut-off
Positive:		OD patient	>	1.2 x OD cut-off	

9 Technical Data

Sample material:	serum
Sample volume:	10 µl of sample diluted 1:101 with 1x sample buffer
Total incubation time:	90 minutes at 20-32°C/68-89.6°F
Calibration range:	0-300 IU/ml
Analytical sensitivity:	3.0 IU/ml
Storage:	at 2-8°C/35-46°F use original vials only.
Number of determinations:	96 tests

10 Performance Data

10.1 Analytical sensitivity

The detection limit of the AESKULISA dsDNA-G has been found as 3 IU/ml.

10.2 Specificity and sensitivity

The microplate is coated with recombinant human dsDNA. No crossreactivities to other autoantigens have been found. For determination of sensitivity and specificity sera of 54 patients suffering from SLE have been tested with the AESKULISA dsDNA-G and an already marketed predicate device. Also 80 healthy donors were tested on the AESKULISA dsDNA-G and were all shown to be negative.

Disease	# of tested sera
SLE	39
MCTD	3
CREST	4
Sjogrens Syndrome	4
Various AI diseases	7

		Predicate device		
		Pos	Non-positive	Total
AESKULISA dsDNA-G	Pos	34	13	47
	Neg	0	7	7
		34	20	54

Clinical specificity: 80/80 = 100%

Clinical sensitivity: 47 / 54 = 87%



10.3 Linearity

Chosen sera have been tested with this kit and found to dilute linearly. However, due to the heterogeneous nature of human autoantibodies there might be samples that do not follow this rule.

Sample No.	Dilution Factor	Measured (IU/ml)	Expected (IU/ml)	Recovery (%)
1	1 / 100	42.9	43.2	99.3
	1 / 200	20.4	21.6	99.4
	1 / 400	9.3	10.8	86.1
	1 / 800	4.9	5.4	90.7
2	1 / 100	179.4	176.0	101.9
	1 / 200	86.4	88.0	98.2
	1 / 400	41.8	44.0	95.0
	1 / 800	19.8	22.0	90.0

10.4 Precision

To determine the precision of the assay, the variability (intra and inter-assay) was assessed by examining its reproducibility on three serum samples selected to represent a range over the standard curve.

Intra-assay		
Sample No.	Mean (IU/ml)	CV (%)
1	199.88	1.00
2	80.14	1.68
3	12.71	5.91

Inter-assay		
Sample No.	Mean (IU/ml)	CV (%)
1	463.3	2.6
2	171.6	2.3
3	13.2	5.9

10.5 Calibration







The AESKULISA dsDNA-G is calibrated against WHO/80. The results are expressed in IU/ml.

11 Literature

Tan EM, Cohen AS, Fries JF, et al. (1982). Revised criteria for the classification of systemic lupus erythematosus. *Arthritis Rheumatism* 25: 1271-1277.

Witte T, Hartung K, Matthias T, Sachse C, Fricke M, Deicher H, Kalden JR, Lakomek HJ, Peter HH, Schmidt RE (1998). Association of IgA anti-dsDNA antibodies with vasculitis and disease activity in systemic lupus erythematosus. *Rheumatol Int* 18: 63-69.

Witte T, Hartung K, Sachse C, Matthias T, Fricke M, Deicher H, Kalden JR, Lakomek HJ, Peter HH, Schmidt RE, SLE study group (1998). IgM anti-dsDNA antibodies in systemic lupus erythematosus: negative association with nephritis. *Rheumatol Int* 18: 85-91.

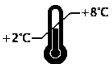IVD	- Diagnosi in vitro	- For in vitro diagnostic use
	- Pour diagnostic in vitro	- Para uso diagnóstico in vitro
	- In Vitro Diagnostikum	- In Vitro Διαγνωστικό μέσο
	- Para uso Diagnóstico in vitro	
REF	° Numero d'ordine	° Catalogue number
	° Référence Catalogue	° Numéro de catálogo
	° Bestellnummer	° Αριθμός παραγγελίας
LOT	° Número de catálogo	
	° Descrizione lotto	° Lot
	° Lot	° Lote
CE	° Chargen Bezeichnung	° Χαρακτηρισμός παρτίδας
	° Lote	
	° Conformità europea	° EC Declaration of Conformity
	° Déclaration CE de Conformité	° Declaración CE de Conformidad
	° Europäische Konformität	° Ευρωπαϊκή συμφωνία
	° Declaração CE de Conformidade	
	° 96 determinazioni	° 96 tests
	° 96 tests	° 96 pruebas
	° 96 Bestimmungen	° 96 προσδιορισμοί
	° 96 Testes	
	° Rispettare le istruzioni per l'uso	° See instructions for use
	° Voir les instructions d'utilisation	° Ver las instrucciones de uso
	° Gebrauchsanweisung beachten	° Λάβετε υπόψη τις οδηγίες χρήσης
	° Ver as instruções de uso	
	° Da utilizzarsi entro	° Use by
	° Utilise avant le	° Utilizar antes de
	° Verwendbar bis	° Χρήση μέχρι
	° Utilizar antes de	
	° Conservare a 2-8°C	° Store at 2-8°C (35-46°F)
	° Conserver à 2-8°C	° Conservar a 2-8°C
	° Lagerung bei 2-8°C	° Φυλάσσεται στους 2-8°C
	° Conservar entre 2-8°C	
	° Prodotto da	° Manufactured by
	° Fabriqué par	° Fabricado por
	° Hergestellt von	° Κατασκευάζεται από
	° Fabricado por	
CO-CAL	° Calibratore cut-off	° Cut off Calibrator
	° Etalon Seuil	° Calibrador de cut-off
	° Grenzwert Kalibrator	° Οριακός ορός Αντιδραστήριο βαθμονόμησης
	° Calibrador de cut-off	
CON+	° Controllo positivo	° Positive Control
	° Contrôle Positif	° Control Positivo
	° Positiv Kontrolle	° Θετικός ορός ελέγχου
	° Controllo positivo	
CON-	° Controllo negativo	° Negative Control
	° Contrôle Négatif	° Control Negativo
	° Negativ Kontrolle	° Αρνητικός ορός ελέγχου
	° Controllo negativo	
CAL	° Calibratore	° Calibrator
	° Etalon	° Calibrador
	° Kalibrator	° Αντιδραστήριο βαθμονόμησης
	° Calibrador	
RC	° Recupero	° Recovery
	° Corrélation	° Recuperado
	° Wiederfindung	° Ανάκτηση
	° Recuperação	
CONJ	° Coniugato	° Conjugate
	° Conjugé	° Conjugado
	° Konjugat	° Σύζευγμα
	° Conjugado	
MP	° Micropiastra rivestita	° Coated microtiter plate
	° Microplaque sensibilisée	° Microplaca sensibilizada
	° Beschichtete Mikrotiterplatte	° Επικαλυμμένη μικροπλάκα
	° Microplaca revestida	
WASHB 50x	° Tampone di lavaggio	° Wash buffer
	° Tampon de Lavage	° Solución de lavado
	° Waschpuffer	° Ρυθμιστικό διάλυμα πλύσης
	° Solução de lavagem	
SUB	° Tampone substrato	° Substrate buffer
	° Substrat	° Tampón sustrato
	° Substratpuffer	° Ρυθμιστικό διάλυμα υποστρώματος
	° Substrato	
STOP	° Reagente bloccante	° Stop solution
	° Solution d'Arrêt	° Solución de parada
	° Stopreagenz	° Αντιδραστήριο διακοπής αντίδρασης
	° Solução de paragem	
SB 5x	° Tampone campione	° Sample buffer
	° Tampon Echantillons	° Tampón Muestras
	° Probenpuffer	° Ρυθμιστικό διάλυμα δειγμάτων
	° Diluente de amostra	