

ORGENTEC Diagnostika GmbH

Carl-Zeiss-Straße 49-51

55129 Mainz - Germany

Phone: +49 (0) 61 31 / 92 58-0

Fax: +49 (0) 61 31 / 92 58-58

Internet: www.orgentec.com



Electronic Instruction For Use: version



917A_3

ORG 917A Anti-Helicobacter pylori IgA

INTENDED PURPOSE

Anti-Helicobacter pylori IgA for Alegria® for Alegria® is an ELISA based test system intended for the quantitative measurement of IgA class antibodies against *Helicobacter pylori* in human serum or plasma samples. This product is intended for professional in vitro diagnostic use only.

The test is used to support the diagnosis of *Helicobacter pylori* infection. *H. pylori* causes a chronic infection of the stomach, which can only be eradicated with a specific antibiotic therapy. Proof of IgG antibodies before treatment onset indicates an existing infection. Some patients develop additional IgA antibodies, which support the IgG findings. After successful treatment with antibiotics, the antibody titres decrease significantly, however, low levels of IgG may persist for up to 12 months.

SYMBOLS USED

| | |
|--|------------------------------------|
| | In vitro diagnostic medical device |
| | Manufacturer |
| | Catalogue number |
| | Sufficient for ... determinations |
| | Batch code |
| | Use by |
| | Temperature limitation |
| | Consult instructions for use |
| | Keep away from sunlight |
| | Do not reuse |
| | Date of manufacture |
| | CE marked according to 98/79/EC |

| | |
|--|----------------------|
| | Alegria® Test Strips |
| | Wash Buffer |
| | System Fluid |
| | Ready to use |

PRINCIPLE OF THE TEST

Reaction wells are coated with a mixture of *H. pylori* whole-cell lysate and recombinant CagA.

The Alegria® assay features barcoded 8-well-microstrips, called Alegria® Test Strips. Each strip is designed for a single determination of one patient sample. The Alegria® Test Strip holds a complete set of reagents. Included are enzyme conjugate, enzyme substrate, sample buffer and a test specific control. Furthermore each strip has two antigen-coated wells which serve as reaction wells for one control and one patient sample.

The determination is based on an indirect enzyme linked immune reaction with the following steps: Antibodies present in positive samples bind to the antigen coated on the surface of the two reaction wells forming an antibody antigen complex. After incubation, a first washing step removes unbound and unspecific bound molecules. Subsequently added enzyme conjugate binds to the immobilized antibody-antigen complex. After incubation, a second washing step removes unbound enzyme conjugate. Addition of enzyme substrate solution results in hydrolysis and color development during incubation. The intensity of the blue color correlates with the concentration of the antibody-antigen-complex and can be measured photometrically at 650 nm.

The Alegria® Test Strip is based on the proprietary SMC®-Technology (Sensotronic Memorized Calibration); information about the assay, analysis and evaluation, and the lot-specific expiry date is contained on the barcode printed on each Alegria® Test Strip.

The Alegria® Test Strip can be used with the diagnostic instrument Alegria® - a fully automated Random Access Analyser. By means of SMC®-Technology data encoded on the barcode are transferred from the Alegria® Test Strip to the instrument and the assay is automatically processed and evaluated. The instrument reads the date of expiry and rejects further processing if the Alegria® Test Strip is out of date.

WARNINGS AND PRECAUTIONS

- All reagents of this kit are intended for professional in vitro diagnostic use only.
- Components containing human serum were tested and found negative for HBsAg, HCV, HIV1 and HIV2 by FDA approved methods. No test can guarantee the absence of HBsAg, HCV, HIV1 or HIV2, and so all human serum based reagents in this kit must be handled as though capable of transmitting infection.
- Bovine serum albumin (BSA) used in components has been tested for BSE and found negative.
- Avoid contact with the substrate TMB (3,3',5,5'-Tetramethyl-benzidine).
- System fluid contains acid, classification is non-hazardous. Avoid contact with skin.
- Control, sample buffer and wash buffer contain sodium azide 0.09% as preservative. This concentration is classified as non-hazardous
- Enzyme conjugate, control and sample buffer contain ProClin 300 0.05% as preservative. This concentration is classified as non-hazardous.


During handling of all reagents, controls and serum samples observe the existing regulations for laboratory safety regulations and good laboratory practice:

- First aid measures: In case of skin contact, immediately wash thoroughly with water and soap. Remove contaminated clothing and shoes and wash before reuse. If system fluid comes into contact with skin, wash thoroughly with water. After contact with the eyes carefully rinse the opened eye with running water for at least 10 minutes. Get medical attention if necessary.
 - Personal precautions, protective equipment and emergency procedures:
Observe laboratory safety regulations. Avoid contact with skin and eyes. Do not swallow. Do not pipette by mouth. Do not eat, drink, smoke or apply makeup in areas where specimens or kit reagents are handled. When spilled, absorb with an inert material and put the spilled material in an appropriate waste disposal.
 - Exposure controls / personal protection: Wear protective gloves of nitril rubber or natural latex. Wear protective glasses. Used according to intended use no dangerous reactions known.
 - Conditions to avoid: Since substrate solution is light-sensitive. Store Alegria® strips in the dark.
 - For disposal of laboratory waste the national or regional legislation has to be observed.
- Observe the guidelines for performing quality control in medical laboratories by assaying controls and/or pooled sera.

CONTENTS OF THE KIT

▽ 24 ORG 917A

ALEGRIA TEST STRIPS 24

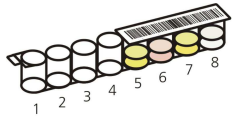

Sufficient for 24 determinations

Alegria® Test Strips are modules of 8 wells each composed of:
Wells 1 + 2: empty and not coated (wells for the sample dilution)

Wells 3 + 4: coated with antigen (reaction wells)

Well 5: Control; yellow; containing test specific antibodies, PBS, BSA, detergent, preservative sodium azide 0.09% and ProClin 300 0.05%.

Well 6: Enzyme Conjugate; light red; containing anti-human IgA antibodies, HRP labelled; PBS, BSA, detergent, preservative ProClin 300 0.05%.

Well 7: Sample Buffer; yellow; containing PBS, BSA, detergent, preservative sodium azide 0.09% and ProClin 300 0.05%.

Well 8: TMB Substrate; clear; containing 3,3', 5,5'- Tetramethylbenzidin.

Code on barcode: **Helico pylori IgA** on printout: **Hpylo-A**

WASH

1x 20 ml

Wash Buffer, containing Tris, detergent, preservative sodium azide 0.09%; 50 x conc.

SYSTEM FLUID

1x 2.5 ml

System Fluid, contains acid; 1000 x concentrate

1

Certificate of Analysis

STORAGE AND STABILITY

- Store test kit at 2-8°C in the dark.
- Do not expose reagents to heat, sun, or strong light during storage and usage.
- Store Alegria® Test Strips sealed and desiccated in the clip bag provided.
- Shelf life of the unopened test kit is 15 months from day of production.
Unopened reagents are stable until expiration of the kit. See labels for individual batch.
- Diluted Wash Buffer and System Fluid are stable for at least 30 days when stored at 2-8°C.
Once transferred to the reagent container we recommend consumption on the same day.

MATERIALS REQUIRED

- Vortex mixer
- Pipettes for 10 µl
- Measuring cylinder for 1000 ml and 2500 ml
- Distilled or deionized water

SPECIMEN COLLECTION, STORAGE AND HANDLING

- Collect whole blood specimens using acceptable medical techniques to avoid hemolysis.
- Allow blood to clot and separate the serum or plasma by centrifugation.
- Test serum should be clear and non-hemolyzed. Contamination by hemolysis or lipemia should be avoided, but does not interfere with this assay.
- Specimens may be refrigerated at 2-8°C for up to five days or stored at -20°C up to six months.
- Avoid repetitive freezing and thawing of serum or plasma samples. This may result in variable loss of antibody activity.
- Testing of heat-inactivated sera is not recommended.

PROCEDURAL NOTES

- Do not use kit components beyond their expiration dates.
- All materials must be at room temperature (20-28°C) prior to use.
- To avoid carryover or contamination, change the pipette tip between samples.

PREPARATION OF REAGENTS

WASH

Dilute the content of the Wash Buffer concentrate (50x) with distilled or deionized water to a final volume of 1000 ml prior to use. Transfer the diluted Wash Buffer into the instrument reagent container. If only one Alegria run is to be performed on one day we recommend transferring only 500 ml diluted Wash Buffer.

SYSTEM FLUID

Dilute the content of the System Fluid concentrate (1000x) with distilled or deionized water to a final volume of 2500 ml prior to use. Transfer the diluted System Fluid into the instrument reagent container.

ALEGRIA TEST STRIPS

Take the required number of Alegria® Test Strips out of the clip bag and let them reach room temperature (20-28°C). Do not remove foil covering the empty wells until you are ready to start the assay.

TEST PROCEDURE

Alegria® Test Strips with SMC® technology are used with the diagnostic instrument Alegria®. Detailed information about operating the instrument can be taken from the Instrument User Manual.


(1) Remove the foil from the empty wells 1 to 4 of the Alegria® Test Strip.

Do not remove foil with printed barcode, covering wells 5 to 8.

(2) Pipette 10 µl undiluted sample at the bottom of well 1.

(3) Insert the strip into the SysTray.

(4) Place loaded SysTrays into the correct position in the Alegria® instrument and start run. All further steps will be done automatically. The test run is completed when the instrument starts printing the results.

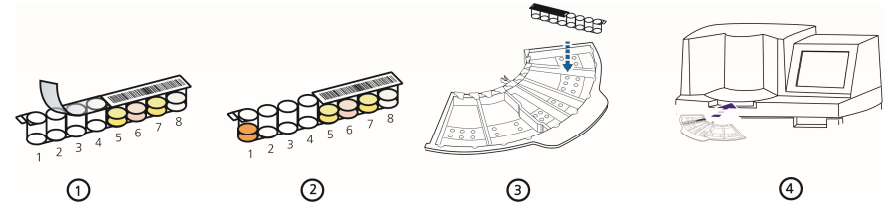

CALIBRATION

This assay system is calibrated in relative arbitrary units, since no international reference preparation is available for this assay.

CALCULATION OF RESULTS

By means of SMC® Technology (Sensotronic Memorized Calibration), all test data are transferred to the system through individual barcodes on the Alegria® Test Strip. Calculation and interpretation of results will be performed automatically.

PERFORMANCE CHARACTERISTICS

Measuring range

The calculation range of this Alegria® assay is: 0 - 200 U/ml

Expected values

The cut-off of this Alegria® assay is: 25 U/ml

Interpretation of results

Negative: < 20 U/ml
Borderline: 20 - 25 U/ml
Positive: > 25 U/ml

LIMITATIONS OF THE PROCEDURE

This assay is a diagnostic aid. A definite clinical diagnosis should not be based on the results of a single test, but should be made by the physician after all clinical and laboratory findings have been evaluated concerning the entire clinical picture of the patient. Also every decision for therapy should be taken individually.

The above reference ranges should be regarded as guidelines only. It is recommended that each laboratory establishes its own normal and pathological ranges for antibodies in patient samples.

A negative result does not rule out an infection, since the serum can be sampled too early for the antibodies to be detectable. A positive result does not rule out the presence of another infectious pathogen as the cause of disease.

Linearity

Three patient samples containing high levels of specific antibody were serially diluted in sample buffer to demonstrate the dynamic range of the assay. Activity for each dilution was calculated by means of SMC® Technology.

| Sample | Dilution | Observed | Expected | O/E |
|--------|----------|----------|----------|-----|
| | | U/ml | U/ml | % |
| 1 | 1:100 | 154.1 | 154.1 | 100 |
| | 1:200 | 78.1 | 77.1 | 101 |
| | 1:400 | 38.0 | 38.5 | 99 |
| 2 | 1:800 | 18.7 | 19.3 | 97 |
| | 1:100 | 114.5 | 114.5 | 100 |
| | 1:200 | 59.4 | 57.3 | 104 |
| 3 | 1:400 | 32.0 | 28.6 | 112 |
| | 1:800 | 16.5 | 14.3 | 115 |
| | 1:100 | 175.4 | 175.4 | 100 |
| . | 1:200 | 82.1 | 87.7 | 94 |
| | 1:400 | 38.8 | 43.9 | 89 |
| | 1:800 | 20.5 | 21.9 | 94 |
| . | 1:1600 | 11.0 | 11.0 | 100 |

Detection limit

The lowest amount of detectable antibody is: 2.7 U/ml

Reproducibility

Intra-assay precision: Coefficient of variation (CV) was calculated for each of three samples from the results of 24 determinations in a single run. Results for precision-within-assay are shown in the table below.

Inter-assay precision: Coefficient of variation (CV) was calculated for each of three samples from the results of 6 determinations in 5 different runs. Results for run-to-run precision are shown in the table below.

| Intra-Assay | | |
|-------------|-------------|------|
| Sample | Mean [U/ml] | % CV |
| 1 | 13.5 | 9.7 |
| 2 | 29.0 | 4.9 |
| 3 | 53.7 | 7.9 |

| Inter-Assay | | |
|-------------|-------------|------|
| Sample | Mean [U/ml] | % CV |
| 1 | 15.2 | 8.5 |
| 2 | 29.1 | 6.8 |
| 3 | 55.6 | 6.8 |

Interfering substances

No interference has been observed with haemolytic (up to 1000 mg/dl) or lipemic (up to 3 g/dl triglycerides) sera or plasma, or bilirubin (up to 40 mg/dl) containing sera or plasma. Nor have any interfering effects been observed with the use of anticoagulants (Citrate, EDTA, Heparine). However for practical reasons it is recommended that grossly hemolyzed or lipemic samples should be avoided.

No interference has been observed in acute EBV infection. Nor have any interfering effects been observed in rheumatic diseases associated with elevated titers of autoantibodies such as rheumatoid factors or antinuclear antibodies.

Seroprevalence

Analysis of 100 healthy blood donors from Germany showed 24 positive results equivalent to 24% seroprevalence.

Study results

| | | Comparative Methods | | |
|------------------------------|--------------|---------------------|-----|-----|
| | | Pos | Neg | |
| Anti-Helicobacter pylori IgA | ORG 917A Pos | 48 | 2 | 332 |
| | Neg | 7 | 275 | |
| | | 55 | 277 | |
| Sensitivity: | | 87.3 | % | |
| Specificity: | | 99.3 | % | |
| Overall agreement: | | 97.3 | % | |

REFERENCES

- Calvet,X., Ramirez Lazaro,M.J., Lehours,P., and Megraud,F., Diagnosis and epidemiology of Helicobacter pylori infection. *Helicobacter*. 2013. 18 Suppl 1: 5-11.
- Fock,K.M., Graham,D.Y., and Malfertheiner,P., Helicobacter pylori research: historical insights and future directions. *Nat Rev Gastroenterol.Hepatol*. 2013. 10: 495-500.
- Gao,L., Michel,A., Weck,M.N., Arndt,V., Pawlita,M., and Brenner,H., Helicobacter pylori infection and gastric cancer risk: evaluation of 15 H. pylori proteins determined by novel multiplex serology. *Cancer Res*. 2009. 69: 6164-6170.
- Leal,Y.A., Flores,L.L., Garcia-Cortes,L.B., Cedillo-Rivera,R., and Torres,J., Antibody-based detection tests for the diagnosis of Helicobacter pylori infection in children: a meta-analysis. *PLoS ONE* 2008. 3: e3751.
- Malfertheiner,P., Megraud,F., O'Morain,C.A., Atherton,J., Axon,A.T., Bazzoli,F., Gensini,G.F., Gisbert,J.P., Graham,D.Y., Rokkas,T., El-Omar,E.M., and Kuipers,E.J., Management of Helicobacter pylori infection--the Maastricht IV/ Florence Consensus Report. *Gut* 2012. 61: 646-664.
- Megraud,F. and Lehours,P., Helicobacter pylori detection and antimicrobial susceptibility testing. *Clin.Microbiol. Rev.* 2007. 20: 280-322.
- Urita,Y., Hike,K., Torii,N., Kikuchi,Y., Kurakata,H., Kanda,E., Sasajima,M., and Miki,K., Comparison of serum IgA and IgG antibodies for detecting Helicobacter pylori infection. *Intern.Med*. 2004. 43: 548-552.
- Wex,T., Venerito,M., Kreutzer,J., Gotze,T., Kandulski,A., and Malfertheiner,P., Serological prevalence of Helicobacter pylori infection in Saxony-Anhalt, Germany, in 2010. *Clin.Vaccine Immunol*. 2011. 18: 2109-2112.


Notice to the user (European Union):

Any serious incident that has occurred in relation to the device shall be reported to the manufacturer and the competent authority of the EU Member State in which the user and/or the patient is established .


Change Control

Former version: *ORG 917A_IFU_EN_QM132288_2014-06-26_2* Reason for revision: *Introduction electronic IFU on homepage*


1

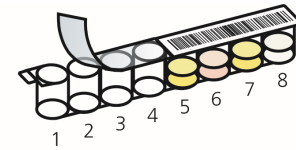

2

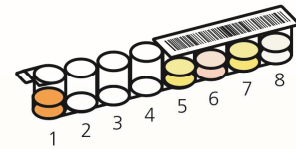

3

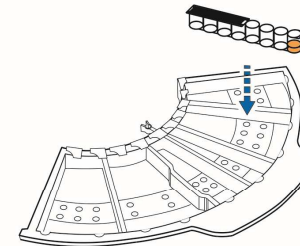

4

