



**DRG<sup>®</sup> Dengue Virus IgM (EIA-3471)**



Revised 3 Oct. 2012 rm (Vers. 16.1)

**RUO** in the USA

**This kit is intended for Research Use Only.**

**Not for use in diagnostic procedures.**

**Intended Use**

The **DRG Dengue Virus IgM Enzyme Immunoassay Kit** provides materials for the **qualitative** and **semiquantitative** measurement of IgM-class antibodies to Dengue Virus in serum.

**PRINCIPLE of the Test**

The **DRG Dengue Virus IgM ELISA Kit** is a solid phase enzyme-linked immunosorbent assay (ELISA)

Specimen samples are diluted with *Sample Diluent* and additionally incubated with *IgG-RF-Sorbent*, containing hyper-immune anti-human IgG-class antibody to eliminate competitive inhibition from specific IgG and to remove rheumatoid factors. This pretreatment avoids false negative or false positive results.

Microtiter wells as a solid phase are coated with Dengue - 2 Virus antigen.

**Diluted sample** specimens and **ready-for-use controls** are pipetted into these wells. During incubation Dengue Virus-specific antibodies of positive specimens and controls are bound to the immobilized antigens.

After a washing step to remove unbound sample and control material horseradish peroxidase conjugated anti-human IgM antibodies are dispensed into the wells. During a second incubation this anti-IgM conjugate binds specifically to IgM antibodies resulting in the formation of enzyme-linked immune complexes.

After a second washing step to remove unbound conjugate the immune complexes formed (in case of positive results) are detected by incubation with TMB substrate and development of a blue color. The blue color turns into yellow by stopping the enzymatic indicator reaction with sulfuric acid.

The intensity of this color is directly proportional to the amount of Dengue Virus-specific IgM antibody in the sample specimen. Absorbance at 450 nm is read using an ELISA microtiter plate reader.

**Warnings and Precautions**

- This kit is for research use only. For professional use only.
- Before starting the assay, read the instructions completely and carefully. Use the valid version of the package insert provided with the kit. Be sure that everything is understood.
- All reagents of this test kit which contain human serum or plasma have been tested and confirmed negative for HIV I/II, HBsAg and HCV by FDA approved procedures. All reagents, however, should be treated as potential biohazards in use and for disposal.
- Avoid contact with Stop Solution containing 0.2 mol/L H<sub>2</sub>SO<sub>4</sub>. It may cause skin irritation and burns.
- TMB substrate has an irritant effect on skin and mucosa. In case of possible contact, wash eyes with an abundant volume of water and skin with soap and abundant water. Wash contaminated objects before reusing them. If inhaled, take the person to open air.
- The microplate contains snap-off strips. Unused wells must be stored at 2 °C to 8 °C in the sealed foil pouch and used in the frame provided
- Pipetting of samples and reagents must be done as quickly as possible and in the same sequence for each step.
- Use reservoirs only for single reagents. This especially applies to the substrate reservoirs. Using a reservoir for dispensing a substrate solution that had previously been used for the conjugate solution may turn solution colored. Do not pour reagents back into vials as reagent contamination may occur.
- Mix the contents of the microplate wells thoroughly to ensure good test results. Do not reuse microwells.
- Do not let wells dry during assay; add reagents immediately after completing the rinsing steps.
- Allow the reagents to reach room temperature (21 °C – 26 °C) before starting the test. Temperature will affect the absorbance readings of the assay. However, values for the specimen samples will not be affected.
- Never pipette by mouth and avoid contact of reagents and specimens with skin and mucous membranes.
- Do not smoke, eat, drink or apply cosmetics in areas where specimens or kit reagents are handled.
- Wear disposable latex gloves when handling specimens and reagents. Microbial contamination of reagents or specimens may give false results.
- Handling should be in accordance with the procedures defined by an appropriate national biohazard safety guideline or regulation.
- Do not use reagents beyond expiry date as shown on the kit labels.
- All indicated volumes have to be performed according to the protocol. Optimal test results are only obtained when using calibrated pipettes and microtiter plate readers.
- Do not mix or use components from kits with different lot numbers. It is advised not to exchange wells of different plates even of the same lot. The kits may have been shipped or stored under different conditions and the binding characteristics of the plates may result slightly different.
- Chemicals and prepared or used reagents have to be treated as hazardous waste according the national biohazard safety guideline or regulation.
- For information on hazardous substances included in the kit please refer to Safety Data Sheets. Safety Data Sheets for this product are available upon request directly from DRG.



Revised 3 Oct. 2012 rm (Vers. 16.1)

RUO	in the USA
-----	------------

### Kit Components

#### Contents of the Kit

1. **Microtiterwells**, 12 x 8 (break apart) strips, 96 wells;  
Wells coated with Dengue - 2 Virus antigen.  
(incl. 1 strip holder and 1 cover foil)
2. **Sample Diluent** \*, 1 vial, 100 mL, ready to use,.  
colored yellow; pH  $7.2 \pm 0.2$ .
3. **IgG-RF-Sorbent** \*, 1 vial, 6.5 mL, ready to use,  
colored yellow;  
Contains anti-human IgG-class antibody.
4. **Pos. Control** \* 1 vial, 2.0 mL, ready to use;  
colored yellow, red cap.
5. **Neg. Control** \*, 1 vial, 2.0 mL, ready to use;  
colored yellow, yellow cap.
6. **Cut-off Control** \*, 1 vial, 2.0 mL, ready to use;  
colored yellow, black cap.
7. **Enzyme Conjugate** \*, 1 vial, 20 mL, ready to use,  
colored red,  
antibody to human IgM conjugated to horseradish peroxidase.
8. **Substrate Solution**, 1 vial, 14 mL, ready to use,  
Tetramethylbenzidine (TMB).
9. **Stop Solution**, 1 vial, 14 mL, ready to use,  
contains 0.2 mol/L H<sub>2</sub>SO<sub>4</sub>,  
Avoid contact with the stop solution. It may cause skin irritations and burns.
10. **Wash Solution** \*, 1 vial, 30 mL (20X concentrated for 600 mL), pH  $6.5 \pm 0.1$   
see „Preparation of Reagents“.

\* contain non-mercury preservative

#### Material required but not provided

- A microtiter plate calibrated reader (450/620nm  $\pm 10$  nm)  
(e.g. the DRG Instruments Microtiter Plate Reader)
- Calibrated variable precision micropipettes
- Incubator 37 °C
- Manual or automatic equipment for rinsing wells
- Vortex tube mixer
- Deionised or (freshly) distilled water



Revised 3 Oct. 2012 rm (Vers. 16.1)

**RUO** in the USA

- Timer
- Absorbent paper

**Storage Conditions**

When stored at 2 °C to 8 °C unopened reagents will retain reactivity until expiration date. Do not use reagents beyond this date.

Opened reagents must be stored at 2 °C to 8 °C. Microtiter wells must be stored at 2 °C to 8 °C. Once the foil bag has been opened, care should be taken to close it tightly again.

Opened kits retain activity for two months if stored as described above.

**Reagent Preparation**

Allow all reagents and required number of strips to reach room temperature prior to use.

**Wash Solution**

Dilute *Wash Solution 1+19* (e.g. 10 mL + 190 mL) with fresh and germ free redistilled water. This diluted wash solution has a pH value of  $7.2 \pm 0.2$ .

Consumption: ~ 5 mL per determination.

Crystals in the solution disappear by warming up to 37 °C in a water bath. Be sure that the crystals are completely dissolved before use.

*The diluted Wash Solution is stable for 4 weeks at 2 °C to 8 °C.*

**Disposal of the Kit**

The disposal of the kit must be made according to the national regulations. Special information for this product is given in the Material Safety Data Sheets.

**Damaged Test Kits**

In case of any severe damage to the test kit or components, DRG has to be informed in writing, at the latest, one week after receiving the kit. Severely damaged single components should not be used for a test run. They have to be stored until a final solution has been found. After this, they should be disposed according to the official regulations.



Revised 3 Oct. 2012 rm (Vers. 16.1)

**RUO** in the USA**SPECIMEN**

Serum can be used in this assay.

Do not use haemolytic, icteric or lipaemic specimens.

*Please note:* Samples containing sodium azide should not be used in the assay.

**Specimen Collection****Serum:**

Collect blood by venipuncture (e.g. Sarstedt Monovette for serum), allow to clot, and separate serum by centrifugation at room temperature. Do not centrifuge before complete clotting has occurred. Samples containing anticoagulant may require increased clotting time..

**Specimen Storage**

Specimens should be capped and may be stored for up to 24 hours at 2 °C to 8 °C prior to assaying.

Specimens held for a longer time should be frozen only once at –20 °C prior to assay. Thawed samples should be inverted several times prior to testing.

**Specimen Dilution**

Prior to assaying each sample specimen is first to be diluted with *Sample Diluent*. For the absorption of rheumatoid factor these prediluted samples then have to be incubated with *IgG-RF-Sorbent*

1. Dilute each sample specimen **1+50** with *Sample Diluent*;  
e.g. 10 µL of specimen + 0.5 mL of *Sample Diluent*. **Mix well.**
2. Dilute this prediluted sample **1+1** with *IgG-RF-Sorbent*  
e.g. 60 µL prediluted sample + 60 µL *IgG-RF-Sorbent*. **Mix well**
3. **Let stand for at least 15 minutes at room temperature, mix well or overnight at 2°C – 8°C and mix well again.**
4. Take 100 µL of these pretreated samples for the ELISA.

**Please note:** Controls are ready for use and must not be diluted!

## Assay Procedure

### General Remarks

- **It is very important to bring all reagents, samples and controls to room temperature before starting the test run!**
- Once the test has been started, all steps should be completed without interruption.
- Use new disposal plastic pipette tips for each standard, control or sample in order to avoid cross contamination
- Absorbance is a function of the incubation time and temperature. Before starting the assay, it is recommended that all reagents are ready, caps removed, all needed wells secured in holder, etc. This will ensure equal elapsed time for each pipetting step without interruption.
- As a general rule the enzymatic reaction is linearly proportional to time and temperature.
- Close reagent vials tightly immediately after use to avoid evaporation and microbial contamination.
- To avoid cross-contamination and falsely elevated results pipette specimen samples and dispense conjugate without splashing accurately to the bottom of wells.

During 37°C incubation cover microtiter strips with foil to avoid evaporation.

### Assay Procedure

Prior to commencing the assay, dilute *Wash Solution*, **prepare specimen samples as described in point 5.3** and establish carefully the **distribution and identification plan** supplied in the kit for all specimens and controls.

Select the required number of microtiter strips or wells and insert them into the holder.

Please allocate at least:

- 1 well (e.g. A1) for the substrate blank,
- 1 well (e.g. B1) for the *Neg. Control*,
- 2 wells (e.g. C1+D1) for the *Cut-off Control*
- 1 well (e.g. E1) for the *Pos. Control*.

It is left to the user to determine controls and specimen samples in duplicate.

### Dispense

- 100 µL** of *Neg. Control* into well B1
- 100 µL** of *Cut-off Control* into wells C1 and D1
- 100 µL** of *Pos. Control* into well E1 and
- 100 µL** of each *pre-treated sample* with new disposable tips into appropriate wells.
- Leave well A1 for substrate blank!



Revised 3 Oct. 2012 rm (Vers. 16.1)

**RUO** in the USA

1. Cover wells with foil supplied in the kit. Incubate for **60 minutes at 37 °C**.
2. Briskly shake out the contents of the wells.  
Rinse the wells **5 times** with diluted *Wash Solution* (**300 µL per well**). Strike the wells sharply on absorbent paper to remove residual droplets.  
**Important note:**  
The sensitivity and precision of this assay is markedly influenced by the correct performance of the washing procedure!
3. Dispense **100 µL Enzyme Conjugate** into each well, **except A1**.
4. Incubate for **30 minutes at room temperature (20 °C to 25 °C)**.  
*Do not expose to direct sun light!*
5. Briskly shake out the contents of the wells.  
Rinse the wells **5 times** with diluted *Wash Solution* (300 µL per well). Strike the wells sharply on absorbent paper to remove residual droplets.
6. Add **100 µL of Substrate Solution** into all wells.
7. Incubate for **exactly 15 minutes at room temperature (20 °C to 25 °C) in the dark**.
8. Stop the enzymatic reaction by adding **100 µL of Stop Solution** to each well.  
Any blue color developed during the incubation turns into yellow.  
**Note:** Highly positive specimen samples can cause dark precipitates of the chromogen!
9. Read the optical density at **450/620 nm** with a microtiter plate reader **within 30 minutes** after adding the *Stop Solution*.

### Measurement

**Adjust** the ELISA microplate or microstrip reader **to zero** using the **substrate blank in well A1**.

If - due to technical reasons - the ELISA reader cannot be adjusted to zero using the substrate blank in well A1, subtract the absorbance value of well A1 from all other absorbance values measured in order to obtain reliable results!

**Measure the absorbance** of all wells **at 450 nm** and record the absorbance values for each control and specimen sample in the distribution and identification plan.

Dual wavelength reading using 620 nm as reference wavelength is recommended.

Where applicable **calculate the mean absorbance values** of all duplicates.

### Results

#### Calculation

##### Mean absorbance value of Cut-off Control [CO]

Calculate the mean absorbance value of the two (2) Cut-off Control determinations (e.g. in C1/D1).

*Example:*  $(0.44 + 0.46) : 2 = 0.45 = CO$



Revised 3 Oct. 2012 rm (Vers. 16.1)

**RUO** in the USA

### Interpretation

**POSITIVE** Sample (mean) absorbance values more than 10 % above CO  
(Mean OD<sub>sample</sub> > 1.1 x CO)

**GREY ZONE** Sample (mean) absorbance values from 10 % above to 10 % below CO  
repeat test 2 - 4 weeks later - with new specimen samples  
(0.9 x CO ≤ Mean OD<sub>sample</sub> ≤ 1.1 x CO)  
Results in the second test again in the grey zone ⇒ **NEGATIVE**

**NEGATIVE** Sample (mean) absorbance values more than 10 % below CO  
(Mean OD<sub>sample</sub> < 0.9 x CO)

### Results in DRG Units [DU]

$$\frac{\text{Sample (mean) absorbance value} \times 10}{\text{CO}} = [\text{DRG Units} = \text{DU}]$$

Example:  $\frac{1.580 \times 10}{0.45} = 35 \text{ DU}$

### Quality Control

It is recommended to use control samples according to state and federal regulations. The use of control samples is advised to assure the day to day validity of results. Use controls at both normal and pathological levels.

It is also recommended to make use of national or international Quality Assessment programs in order to ensure the accuracy of the results.

In this case, please check the following technical areas: Pipetting and timing devices; photometer, expiration dates of reagents, storage and incubation conditions, aspiration and washing methods.

After checking the above mentioned items without finding any error contact your distributor or DRG directly.

### Legal Aspects

#### Reliability of Results

The test must be performed exactly as per the manufacturer's instructions for use. Moreover the user must strictly adhere to the rules of GLP (Good Laboratory Practice) or other applicable national standards and/or laws. This is especially relevant for the use of control reagents. It is important to always include, within the test procedure, a sufficient number of controls for validating the accuracy and precision of the test.



# DRG



**DRG<sup>®</sup> Dengue Virus IgM (EIA-3471)**



**Revised 3 Oct. 2012 rm (Vers. 16.1)**

**RUO in the USA**

### **Liability**





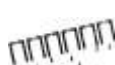




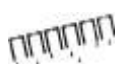

Any modification of the test kit and/or exchange or mixture of any components of different lots from one test kit to another could negatively affect the intended results and validity of the overall test. Such modification and/or exchanges invalidate any claim for replacement.

In the event of any claim, the manufacturer's liability is not to exceed the value of the test kit. Any damage caused to the test kit during transportation is not subject to the liability of the manufacturer.

### **REFERENCES**

1. Hayes E.B., Gubler D.J., Dengue and dengue hemorrhagic fever, *Pediatric Inf.Dis.J.* 11:311-317, 1992
2. CDC: Imported Dengue-United States 1993-1994, *J.Am.Med.Assn.* 274:113, 1995
3. Rigau-Perez J.G., Gubler D.J., Vorndam A.V. et al, Dengue surveillance-United States 1986-1992, *Morbidity and Mortality Weekly Report* 43:7-19, 1994

### Short Instructions for Use






	<p>All reagents and specimens must be allowed to come to room temperature (18-25°C) before use.</p>
	<p>Leave well A1 for substrate Blank. Dispense 100 µl of Controls into appropriate wells.</p>
	<p>Dispense 100 µl of sample into selected wells. <b>(Please note special sample treatment, point 5.3!)</b></p>
	<p>Cover wells with foil. Incubate for <b>60 minutes</b> at 37 °C.</p>
	<p>Briskly shake out the contents of the wells.</p>
	<p>Rinse the wells <b>5 times</b> with diluted Wash Solution (300 µl per well).</p>
	<p>Strike the wells sharply on absorbent paper to remove residual droplets.</p>
	<p>Dispense 100 µl of Enzyme-Conjugate into each well.</p>
	<p>Incubate for <b>30 minutes</b> at room temperature.</p>
	<p>Briskly shake out the contents of the wells.</p>
	<p>Rinse the wells <b>5 times</b> with diluted Wash Solution (300 µl per well).</p>

## DRG<sup>®</sup> Dengue Virus IgM (EIA-3471)



Revised 3 Oct. 2012 rm (Vers. 16.1)

**RUO** in the USA

	<p>Strike the wells sharply on absorbent paper to remove residual droplets.</p>
	<p>Add 100 µl of Substrate Solution to each well.</p>
 <p>15 min</p>	<p>Incubate for <b>15 minutes</b> at room temperature.</p>
	<p>Stop the reaction by adding 100 µl of Stop Solution to each well.</p>
	<p>Determine the absorbance of each well at 450 nm.</p>

Rev. 9/28/12cc